


# Ovens and Incubators

## Incubators



Mini

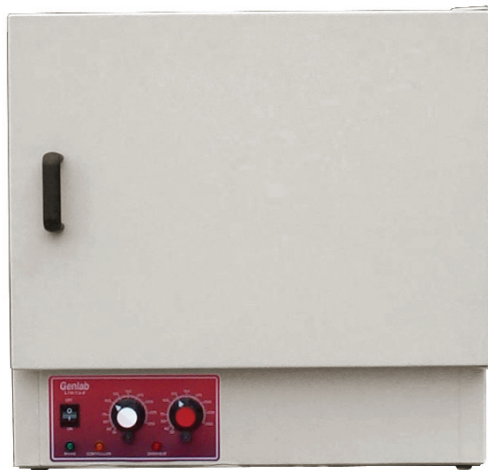

Model INC/100

A range of good quality general purpose horizontal and vertical style incubators with a range of options to suit the needs of every application. The exterior and doors are of easy to clean coated steel. Internally, the chamber is of mild steel, coated with aluminium, and fitted with a thermometer clip. Removable shelf grids are of chrome-plated wire. A direct-reading thermostat for control plus independent overheat thermostat for safety are fitted as standard, with heat and overheat indicators. Both thermostats have scales calibrated for accuracy and tamper-proof locks.

The incubators have a temperature range from ambient +5°C (without fan, or +8°C with fan) to 100°C (with steel door: or 60°C with transparent door). All models require a power supply of 220/240 volts, 50Hz, single phase and feature a power isolation switch with indicator.

### Vertical Models:

Product No	Model	Capacity-Lt	Internal H x W x D-mm	External H x W x D-mm*	Shelves	No.positions	Nett weight-kg
11239/01	MINI/6	6	150 x 230 x 180	390 x 320 x 320	1	1	7
11239/02	MINI/18	18	260 x 260 x 270	500 x 360 x 400	2	2	14
11239/03	MINI/30	30	240 x 360 x 350	490 x 450 x 470	2	2	20
11239/04	MINI/40	40	320 x 360 x 350	570 x 450 x 470	2	3	23
11239/05	MINI/50	40	330 x 490 x 330	580 x 610 x 450	2	3	26
11239/06	MINI/75	75	330 x 490 x 460	580 x 610 x 580	2	3	34
11239/07	MINI/100	100	450 x 490 x 460	700 x 610 x 580	3	4	44
11239/08	MINI/125	125	550 x 490 x 460	800 x 610 x 580	3	5	52
11239/09	MINI/150	150	550 x 490 x 540	800 x 610 x 660	3	5	60
11239/10	MINI/175	175	670 x 490 x 540	920 x 610 x 660	4	6	68
11239/11	MINI/200	200	750 x 490 x 540	1000 x 610 x 660	4	7	76

\*NB above external depths are for models WITHOUT fans. Where a fan is fitted, the depth will increase by 120mm

Product No	Model	Capacity-Lt	Internal H x W x D-mm	External H x W x D-mm*	Shelves	No.positions	Nett weight-kg
11240/01	INC/75	75	420 x 530 x 350	550 x 850 x 460	2	3	39
11240/02	INC/100	100	420 x 530 x 460	550 x 850 x 570	3	4	48
11240/03	INC/125	125	520 x 530 x 460	650 x 850 x 570	3	5	57
11240/04	INC/150	150	520 x 530 x 540	650 x 850 x 650	3	5	66
11240/05	INC/175	175	640 x 530 x 540	770 x 850 x 650	4	6	75
11240/06	INC/200	200	720 x 530 x 540	850 x 850 x 650	4	7	85

\*NB above external depths are for models WITHOUT fans. Where a fan is fitted, the depth will increase by 120mm